



Omnia Stapling Platform

Ultraflex 90

90° Articulation Reload



Omnia
Platform



90°
Articulation



Pro Height
Staples



30 / 45 / 60 mm
Staple Length

For LATAM Use Only
Available products subject to
local regulations and registrations

Product Features



Industry-first 90° Articulation

Witness a monumental pivot in stapling technology.

Offering unparalleled 90-degree articulation, Ultraflex 90 reloads provide surgeons with exceptional freedom of maneuverability, allowing them to reach farther and access difficult areas with unprecedented precision.


Specifications

Product Code	Staple Line Length	Trocar Compatibility (mm)	Color	Staple Height (Open/Closed) (mm)	Tissue Thickness	Number of staple rows	QTY/ BOX
IDT30TAN-R IDT30TAN-L	30mm	12	Tan	2.0/2.5/3.0 0.75/1.0/1.25	Vascular/ Thin	6	6
IDT45TAN-R IDT45TAN-L	45mm						
IDT60TAN-R IDT60TAN-L	60mm						
IDT30PUL-R IDT30PUL-L	30mm		Purple	3.0/3.5/4.0 1.25/1.5/1.75	Medium/ Thick		
IDT45PUL-R IDT45PUL-L	45mm						
IDT60PUL-R IDT60PUL-L	60mm						
IDT45BLK-R IDT45BLK-L	45mm	15	Black	4.0/4.5/5.0 1.75/2.0/2.2	Extra Thick		
IDT60BLK-R IDT60BLK-L	60mm						


Optimized for Lower Anterior Resection (LAR)

The Ultraflex 90 excels in confined, narrow spaces, making it particularly suited for Lower Anterior Resection (LAR) procedures. Surgeons can maintain a safe distance from the tumor site while minimizing the risk of damaging surrounding vital structures.

Lower Reach in Confined Spaces

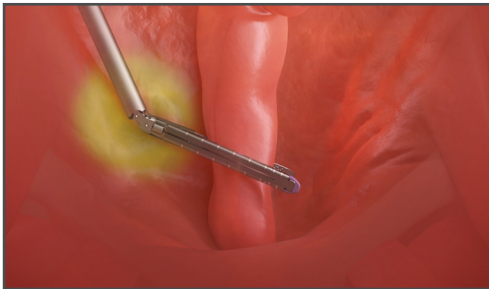

VS

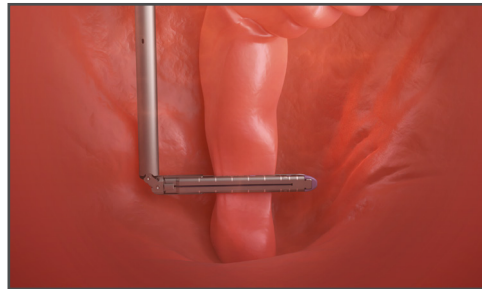


Conventional Reloads

Conventional reloads with 45-degree articulation can be limiting when positioning at lower angles.


Ultraflex 90 Reloads

Ultraflex 90's unique 90-degree articulation allows surgeons to reach deeper into the pelvic region, ensuring a clean transection away from the tumor site with improved accessibility.

Reduced Leakage Risk

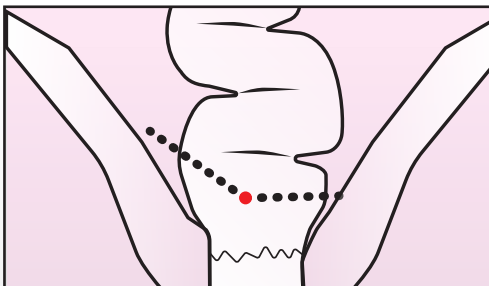

VS

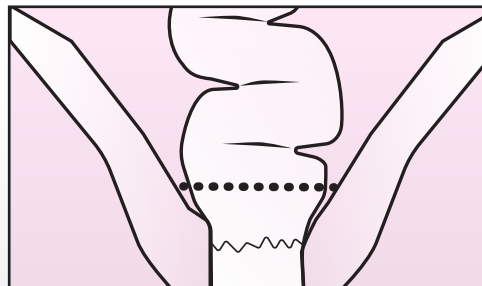

Conventional Reloads

Conventional reloads often require multiple firings, resulting in crossing staple lines that significantly increase the risk of leakage.

Ultraflex 90 Reloads

Ultraflex 90 achieves a clean transection with a single firing, reducing the risk of leakage caused by crossed staple lines.

Other Applications

The Ultraflex 90 offers surgeons increased flexibility across a wide range of surgeries. Tissue approximation is made easier, with less tension and more access than ever before.



Colorectal



Thoracic

Contacts

International Headquarters

16 Science Park Drive
#04-03 DNV GL
Technology Centre
Singapore 118227

T +65 6865 9879

Importer



Distributor

